


LeviKam Standard Sling

LST-XS/S/M/L/XL-1



The LeviKam Standard sling is compatible with hoists for people with disabilities. It features high-quality workmanship and a well-thought-out design. It is made of durable, safe polyester mesh that is highly water-permeable and dries quickly. The LeviKam Standard sling allows the user to be transported in a seated position while providing back support. The user's position can be adjusted to suit their needs thanks to the four-point, adjustable suspension. A sewn-in maneuvering handle allows the user to be rotated while hanging.

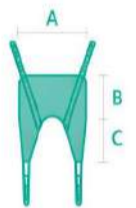

LeviKam Standard	XS size	S size	M size	L size	XL size
A size	66	78	90	101	112
B size	50	51	52	54	56
C size	51	55	59	76	88

*The values given are approximate due to differences in body structure.

Main advantages of slings:

- Maximum load capacity: 200 kg,
- Optional seats with head support,
- Multi-level adjustment of the vest attachment points,
- Wide range of sizes available: from XS to XL,
- Slings feature a physiological cutout to facilitate toileting and hygiene procedures,
- Seats are available in two versions: a standard bath seat (made of polyester mesh) or a special comfort seat (made of solid polyester material, sewn upon request).

Stiffened head support L-UP-1

A practical and comfortable head stabilizing insert, designed to complement LeviKam vests, which do not normally have head support.



Instructional videos

Putting on slings in a sitting position:



Putting on slings in a lying position:

